

# Cool Cube™ 03 at Fridge Temps (Blood Products)

Model: BT-03




The Cool Cube™ 03 at Fridge Temperatures for Blood Products is not an ordinary transport cooler. It is an easy-to-carry, qualified container and pack-out system that keeps temperature sensitive product (whole blood, RBCs, liquid plasma) chilled between 1° and 10°C for over 2.5 days – without using ice or electricity! This cooler is ideal for transport, but also excels at protecting valuable product during off-site use and power outages, regardless of hot or freezing conditions. The simple, validated, patented technology is the best compact thermal protection product available to maintain the cold chain!




## Validated Pack-Out System


Remove PCM Panels

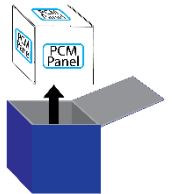


Prep PCM Panels

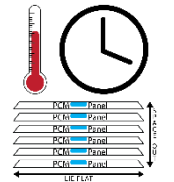

Assemble PCM Panels

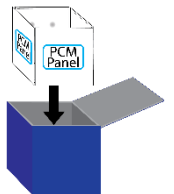

Pack-Out Product

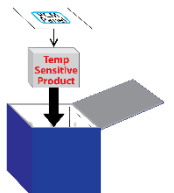

Close Bag



### Features

- 6 patented PCM panels manage the temp
- Temp-Shield™ vacuum insulated panel (VIP) system
- Rugged, blue carrying case with a shoulder & hand strap
- Validated under ISTA Test Procedure 7D (Int'l Safe Transit Assoc)

### Specifications

- Temperature Range*  
1° to 10°C (34° to 50°F)
- Performance*  
70+ Hours
- Payload Area*  
5¾" x 5¾" x 5¾"
- Exterior Dimensions & Weight*  
11" x 11" x 11" @ 11lbs

- 1 Adjustable shoulder strap with pad is removable.
- 2 Double zipper closure for easy access.
- 3 Clear external pockets for documents and/or Digital Data Logger (DDL).
- 4 Six removable, temperature-controlling, PCM panels ensure a simple, quick & effective pack-out.

Go to Web Page >>



info@VeriCorMed.com  
(866) 469-6019

For more information & pricing: [www.VeriCorMed.com](http://www.VeriCorMed.com)

