

Cool Cube™ 03 at Room Temps (FFPE, platelets, biospecimens)

Model: RT-03




The Cool Cube™ 03 at Room Temperatures is not an ordinary transport cooler. It is an easy-to-carry, qualified container and pack-out system that keeps temperature sensitive product (FFPE, platelets, biospecimens) between 15° and 25°C for over 2.5 days – without using ice or electricity! This cooler is ideal for transport, but also excels at protecting valuable product during off-site use and power outages, regardless of hot or freezing conditions. The simple, validated, patented technology is the best compact thermal protection product available to maintain the cold chain!




Validated Pack-Out System


Remove PCM Panels

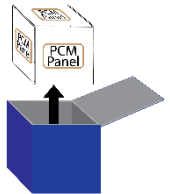


Prep PCM Panels

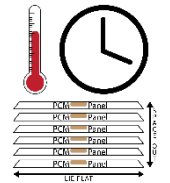

Assemble PCM Panels

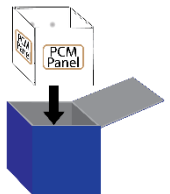

Pack-Out Product

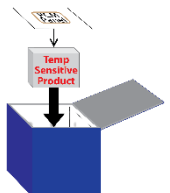

Close Bag



Features

- 6 patented PCM panels manage the temp
- Temp-Shield™ vacuum insulated panel (VIP) system
- Rugged, blue carrying case with a shoulder & hand strap
- Validated under ISTA Test Procedure 7D (Int'l Safe Transit Assoc)

Specifications

- Temperature Range*
15° to 25°C (59° to 77°F)
- Performance*
91+ Hours (summer)
- Payload Area*
5¾" x 5¾" x 5¾"
- Exterior Dimensions & Weight*
11" x 11" x 11" @ 11lbs

- 1 Adjustable shoulder strap with pad is removable.
- 2 Double zipper closure for easy access.
- 3 Clear external pockets for documents and/or Digital Data Logger (DDL).
- 4 Six removable, temperature-controlling, PCM panels ensure a simple, quick & effective pack-out.

Go to Web Page >>



info@VeriCorMed.com
(866) 469-6019

For more information & pricing: www.VeriCorMed.com

