

Cool Cube™ 08 at Refrigerator Temps (Vaccine & Medicine)

Model: VT-08







The Cool Cube™ 08 at Refrigerator Temperatures for Vaccine & Medicine is not an ordinary transport cooler. It is an easy-to-carry, qualified container and pack-out system that keeps temperature sensitive product (fresh vaccine, medicine) chilled between 2° and 8°C for over 3 days – without using ice or electricity! This cooler is ideal for transport, but also excels at protecting valuable product during off-site use and power outages, regardless of hot or freezing conditions. The simplified, validated, patented technology is the best compact thermal protection product available to maintain the cold chain!



4 Crēdo™ inside

- 1 Adjustable shoulder strap with pad is removable.
- 2 Double zipper closure for easy access.
- 3 Clear external pockets for documents and/or Digital Data Logger (DDL).
- 4 Six removable, temperature-controlling, PCM panels ensure a simple, quick & effective pack-out.

Qualified Pack-Out System

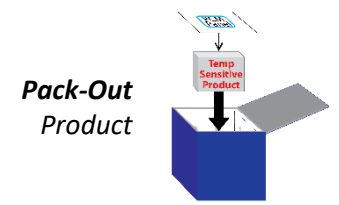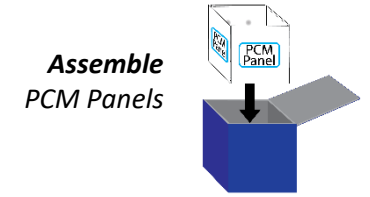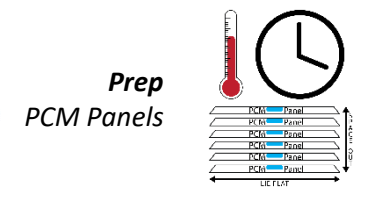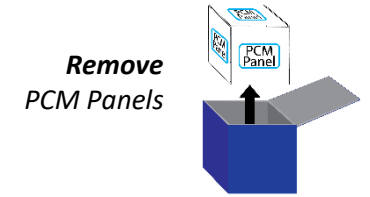

Features

- 6 patented PCM panels manage the temp
- Temp-Shield™ vacuum insulated panel (VIP) system
- Rugged, blue carrying case with a shoulder & hand strap
- Validated under ISTA Test Procedure 7D (Int'l Safe Transit Assoc)

Specifications

Temperature Range
2° to 8°C (36° to 46°F)

Performance
76+ Hours

Payload Area
11¼" x 6½" x 6½"

Exterior Dimensions & Weight
16" x 11½" x 11½" @ 16lbs

[Go to Web Page >>](#)



info@VeriCorMed.com
(866) 469-6019

For more information & pricing: www.VeriCorMed.com

