

Cool Cube™ 96 at Fridge Temps (Blood Products)

Model: BT-96







The Cool Cube™ 96 at Fridge Temperatures for Blood Products is not an ordinary transport cooler. It is a lockable, wheeled, validated container and pack-out system that keeps temperature sensitive product (whole blood, RBCs, liquid plasma) chilled between 1 and 10 °C for over 5 days (or between 1 and 6 °C for over 4 days) – without using ice or electricity! This cooler is ideal for transport, but also excels at protecting valuable product during off-site use and power outages, regardless of hot or freezing conditions. The simple, validated, patented technology is the best compact thermal protection product available to maintain the cold chain!



- 1 Folding handle for ease of movement.
- 2 Holder and flexible grommet enables the use of an external temp monitor.
- 3 Multi-directional locking wheel system.
- 4 Six removable, temperature-controlling, PCM panels ensure a simple, quick & effective pack-out.

Validated Pack-Out System

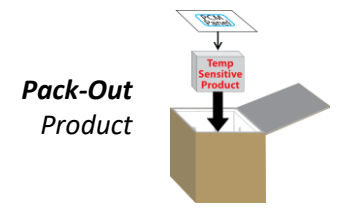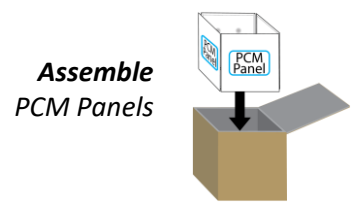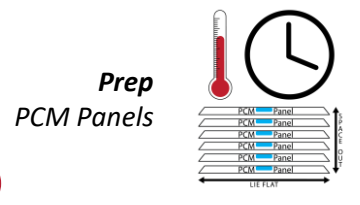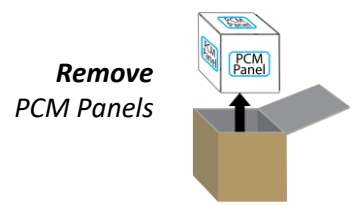

Features

- 6 patented refrigerator PCM panels manage the temp
- Temp-Shield™ vacuum insulated panel (VIP) system
- Rugged, military spec case with large lifting handles & wheels
- Validated under ISTA Test Procedure 7D (Int'l Safe Transit Assoc.)

Specifications

Temperature Range / Duration
 1 to 6 °C (34 to 42 °F) / 112+ hrs.
 1 to 10 °C (34 to 50 °F) / 128+ hrs.

Payload Area
 18" x 18" x 18"

Exterior Dimensions & Weight
 26½" x 26½" x 28¼" @ 111 lbs.

Go to Web Page >>



info@VeriCorMed.com
 (608) 526-6901

For more information & pricing: www.VeriCorMed.com

