

Traceable® Fridge/Freezer Temperature Monitoring Kit

Model: TM-5650

The Traceable® Fridge/Freezer Temperature Monitoring Kit reliably monitors the temperature and alerts you of excursions without opening the door/container. This calibrated digital thermometer has a detachable, GRAS, glycol-encased probe (recommended by the CDC) to insulate the readings from rapid temperature changes (i.e. when a door is opened). The integrated alarm bar lights up and buzzes if an excursion occurs (the alarm light also glows orange within 1 degree of alarm settings). The triple display simultaneously shows minimum, maximum, and current temperatures for effective monitoring. VeriCor is a proud distributor of Traceable® products.



Features

- Traceable to NIST for accuracy
 - Certificate included
- Records high/low readings
- Signals out-of-range conditions (visual & audible)
- 10 foot micro-thin cable
- Integrated stand & magnet

Specifications

Temperature Range
-50° to 70°C (-58° to 158°F)

Resolution
1°


Accuracy
± 1° C

Dimensions & Weight
2¾" x 5¼" x ¾" @ 7oz

Front Panel Open



Back Features

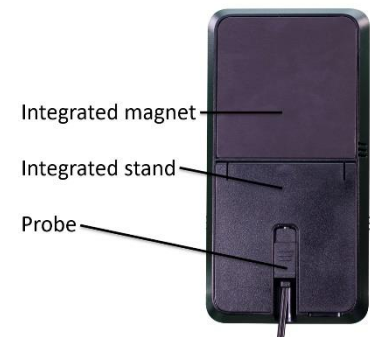

Go to Web Page >>



info@VeriCorMed.com
(608) 526-6901

For more information & pricing: www.VeriCorMed.com

VERICOR LLC
Response-Friendly®